

# Lisa<sup>®</sup>

A sterile state is difficult to maintain. It is a condition that lasts only as long as an absolutely resistant and hermetically sealed protective system protects the sterile object from the environment.

DISTILLER

TESTER

## Mseal 12


PRINTER

STERILIZER

SEALER

# W&H

The **Lisa<sup>®</sup>** Mseal 12 sealer produces perfect seals guaranteeing that bags are hermetically sealed.

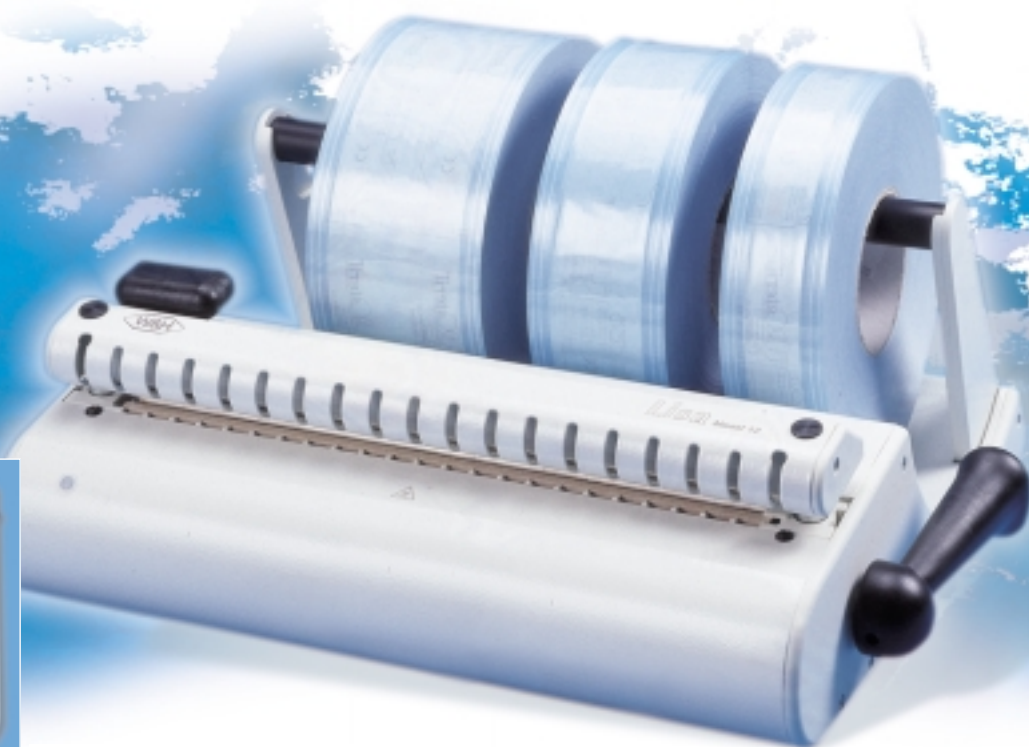

Wall installation.

### Savings:

- The sealer can take various widths of roles which makes it possible to suit the size of the bag to that of the object to be sterilized.
- Low power consumption.

### Effectiveness:

- The temperature of the heating element is controlled electronically thus guaranteeing unbroken seals 12 mm wide.
- 2 to 3 seconds are enough to seal the sachet.
- The end of the operation is indicated by an acoustic signal.

### Safety:

- The protective housing prevents any contact with the heating element and the moving cutter.
- Perfect quality seals with no risk of burning the bags.

### Technical data:



Power consumption	100 W / 0,4 A
Size	L: 480 mm / H: 198 mm / D: 385 mm
Weight	4 kg
Type of seal	constant temperature
Width of seal	12 mm
Max bag width	310 mm
Sealing process time	2 to 3 seconds

1 YEAR  
WARRANTY